

## LABTECH INTERNATIONAL


<b>Item Name :</b>	Demonstration Transformer
<b>Item Code :</b>	13720
<b>Item Description :</b>	

### Product Details:

Phase	Single Phase
Voltage Rating	Input Voltage: 220 Volt and Output Voltage: 12 Volts
Cooling Type	Dry Type/Air Cooled
Power(VA)	12 Volt

Fitted on board with one main coil of 220 volt and two interchangeable secondary coil giving 6 and 12 volts AC 2Amps. Assembly and interchanging of coil can be done easily. input and output leads brought out on terminal.

**Item Image :**

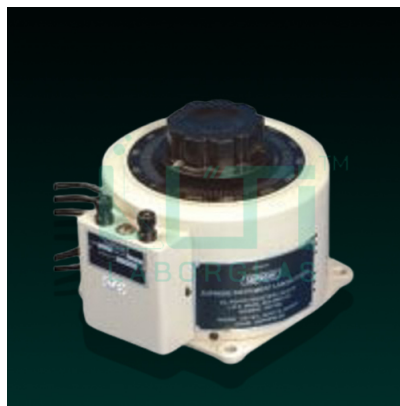

11-A, Herman City Road, Khudda Khurd, Ambala Cantt, Haryana 133006

+91-9802013010, +91-9896317240

info@labtechexport.com, <http://www.labtechexport.com>